Instructions for use

MK-dent High-Speed Turbine

Version: 2026-09-01-RevC

Language Index

English Name	Native Name	ISO 639-1	Page
English	English	en	3
French	Français	fr	34
Spanish	Español	es	66
German	Deutsch	de	97
Italian	Italiano	it	129
Chinese	中文	zh	160
Swedish	Svenska	sv	192
Hungarian	Magyar	hu	223
Slovak	Slovenčina	sk	254
Czech	Česky	cs	285
Turkish	Türkçe	tr	316
Norwegian	Norsk	no	347
Bulgarian	Български	bg	378
Lithuanian	Lietuvių	lt	409
Portuguese (Portugal)	Português	pt	440
Arabic	العربية	ar	472
Greek	Αγγλικά	el	503

Table of Contents

1	About these instructions for use	6
1.1	Explanation of warnings	6
1.2	Compliance	6
2	Intended purpose	7
2.1	Intended use	7
2.2	Indications	7
2.3	Contraindications	7
2.4	Intended patient groups	7
2.5	Intended users	7
2.6	User environment	7
3	Safety and requirements	8
3.1	General safety regulations	8
3.2	Reporting obligation	8
3.3	Side effects	8
3.4	Risk patient groups	8
4	Product description	9
4.1	Description	9
4.2	Combination with other products	10
4.3	Product components	10
4.4	Scope of delivery	10
4.5	Required materials that are not included	11
4.6	Configurations	11
5	Accessories, consumables and spare parts	12
5.1	Accessories	12
5.2	Consumables	13
5.3	Spare parts	13
6	Installation	14
6.1	Requirements for the environment	14
6.2	Unpacking	14
6.3	Assembly	14
6.4	Preparation for use	14
6.4.1	Check and setting of the compressed air	15
6.4.2	Connecting to the coupling system	15
6.4.2.1	Connecting the instrument – turbine with direct connection	15
6.4.2.2	Connecting the instrument – turbine with coupling connection	16

6.5	Startup and operation	16
6.5.1	Before use	16
6.5.2	Checking for correct operation	16
6.5.3	Testing the water quantity	16
6.5.4	Using of rotating instruments	17
7	Treatment	18
7.1	Patient, operator and assistant protection	18
7.2	Removing rotating instruments	18
8	Processing according to ISO 17664-1	19
8.1	Pre-treatment at the point of use	20
8.1.1	Preparation for cleaning	20
8.1.2	Manual pre-cleaning of the spray nozzle	21
8.1.3	Manual pre-cleaning of the instrument	21
8.1.4	Manual pre-cleaning of the glass rod of turbines with light	22
8.2	Manual processing	23
8.2.1	Manual cleaning	23
8.2.2	Manual disinfection, rinsing and drying	23
8.3	Automatic processing	24
8.3.1	Automatic cleaning, disinfection and rinsing	24
8.3.2	Automatic drying	24
8.3.3	Maintenance, checks and testing	25
8.4	Sterilization in an autoclave	25
8.4.1	Packaging	25
8.4.2	Sterilization in a steam sterilizer according to EN 13060 / EN ISO 17665	26
8.4.3	Sterilization parameters	26
8.4.4	After sterilization	26
8.5	Storage	27
9	Maintenance	27
9.1	Regular maintenance and checks	27
9.2	Cleaning in case of problems	27
10	Transport and Disposal	28
11	Troubleshooting	28
12	Repair	29
13	Technical data	30
14	Explanation of symbols	31
15	Guarantee and warranty	32
15.1	Guarantee	32

15.2 Legal warranty 32

15.3 Conditions..... 33

15.4 Utilization 33

15.5 Exclusion of liability 33

1 About these instructions for use

These instructions for use apply to all MK-dent turbines (see section **4 Product description**).

All products listed here are only intended for dental treatment in the field of dentistry. Any kind of misuse or modification of the product is not permitted and may lead to hazards. All products listed here may only be used by qualified users in accordance with **2.5 Intended users** for the intended purpose and in compliance with the instructions for use.

1.1 Explanation of warnings

These instructions for use are intended to describe the safe and effective operation of the product to the user. Before operating the product, you must read these instructions for use, observe all safety instructions and warnings and strictly adhere to them. Pay particular attention to all information and procedures described in section **3 Safety and requirements**.

These instructions for use are part of the product and must be kept within reach of the product so that they are always accessible.

To avoid personal injury and damage to property, observe the warnings and safety information contained in this document. These are specially marked:

WARNING

A WARNING alerts you to a potentially serious outcome, unwanted result or safety risk. Failure to observe a warning can result in the death or serious injury of the user or patient.

CAUTION

CAUTION alerts you to situations requiring particular care for the safe and effective use of the product. Disregard for a precautionary measure can result in minor or moderately severe personal injuries or damages to the product or other property and may potentially lead to a remote risk of serious injuries and/or environmental pollution.

NOTE

NOTE

This identifies special tips, such as help to the user or for the improvement of a work process.

1.2 Compliance

If you have any further questions regarding the applicable national or international standards, please contact the manufacturer:

MK-dent GmbH, Marie-Curie-Str. 2, 22941 Bargteheide,
Germany
+49 4532 40049-0
info@mk-dent.com

2 Intended purpose

2.1 Intended use

High-speed handpieces (**turbines**, max. 250,000 – max. 420,000 min⁻¹): drives for rotating instruments for the removal of carious material, for the preparation of cavities and crowns and for the removal of fillings.

2.2 Indications

Diseases of the periodontium, teeth and jaws (ICD-10 codes).

Turbines are drives for rotating instruments used for the removal of caries, the preparation of cavities and crowns and the removal of fillings.

2.3 Contraindications

There are currently no known contraindications.

2.4 Intended patient groups

The decision whether to use a turbine on a patient is made by the treating dentist. Depending on the assessment of the user, the products can be used on all patients who require dental treatment.

2.5 Intended users

MK-dent turbines are intended exclusively for use by accredited dentists. Secondary users are the practice staff who are responsible for transport, processing and cleaning.

2.6 User environment

Application takes place in dental practices or dental facilities. Users wear protective gloves. The peripherals used fulfil the requirements of ISO 7494.

3 Safety and requirements

3.1 General safety regulations

WARNING

- This product is supplied non-sterile. Prepare it before first use in accordance with section **8 Processing according to ISO 17664-1**.
- For proper reprocessing, reprocess the device after each patient in accordance with section **8 Processing according to ISO 17664-1**. Safe and effective reprocessing after multiple use without reprocessing the device leads to mixed infections and prevents proper processing due to an atypical accumulation of contamination, closures and cleaning agent residues – especially in cavities and supply lines.
- This product may be used only by qualified personnel and only for dental treatment.
- Never disassemble or modify this product.
- This instrument must not be repaired by the user in the event of a malfunction. Contact your dealer or the MK-dent repair.
- Please do not use the product if the performance of the product deteriorates or in case of malfunctions.
- Do not use this product if it has not been serviced properly.
- Observe the precautionary measures for use in order to ensure correct and safe functioning.

CAUTION

- The product may only be used in the field of dentistry with a dental treatment unit that has a CE marking or is otherwise certified/approved according to the valid national legal regulations and that is provided with an appropriate water and compressed air supply in compliance with the applicable standards.
- Test for proper functioning (collet holding force, spray water flow rate and heat development) before each use and inspect the instrument visually for external damage and contamination.
- If the amount of spray water is too low, the instrument may overheat and cause injuries.

3.2 Reporting obligation

In the event of serious deterioration in the patient's health during the use of this product, please report this to your specialized dealer or to us as the manufacturer and to the national competent authority.

Our contact details are included in these instructions for use. The contact details of the notified body will also be provided to you upon request.

3.3 Side effects

There are currently no known side effects.

3.4 Risk patient groups

There are currently no known risk patient groups.

4 Product description

4.1 Description

MK-dent turbines are available in the following variants: PRIME LINE, CLASSIC LINE, ECO LINE, BASIC LINE

The compact shape of the turbine makes switching between treatment units easy and requires very little space for storage.

Thanks to the quick-coupling system, swapping turbines is easy, enabling immediate use. No complicated installation is required, allowing for simple integration in your treatment system.

Driven by compressed air from the treatment unit, the turbine has no effect on other medical devices, nor is it affected by electromagnetic waves.

Article no.	Light	Standard head	Small head	Mini head	Direct connection	Quick-coupling system
HB12	-	X	-	-	4-hole	-
HB14	-	-	-	X	4-hole	-
HB15L	X	X	-	-	5-hole	-
HB16L	X	-	-	X	5-hole	-
HB21K	-	X	-	-	-	KaVo
HB21KL	X	X	-	-	-	KaVo
HB23K	-	-	-	X	-	KaVo
HB23KL	X	-	-	X	-	KaVo
HE11	-	X	-	-	2-hole	-
HE12	-	X	-	-	4-hole	-
HE14	-	-	-	X	4-hole	-
HE15L	X	X	-	-	5-hole	-
HE16L	X	-	-	X	5-hole	-
HE17	-	X	-	-	4-hole	-
HE17B	-	X	-	-	2-hole	-
HE18	-	-	X	-	4-hole	-
HE20K	-	X	-	-	-	KaVo
HE20KL	X	X	-	-	-	KaVo
HE21K	-	X	-	-	-	KaVo
HE21KL	X	X	-	-	-	KaVo
HE21N	-	X	-	-	-	NSK
HE21NL	X	X	-	-	-	NSK
HE21SL	X	X	-	-	-	Sirona
HE21W	-	X	-	-	-	W&H
HE21WL	X	X	-	-	-	W&H
HE22K	-	-	X	-	-	KaVo
HE22KL	X	-	X	-	-	KaVo
HE22N	-	-	X	-	-	NSK
HE22NL	X	-	X	-	-	NSK
HE22SL	X	-	X	-	-	Sirona
HE22W	-	-	X	-	-	W&H
HE22WL	X	-	X	-	-	W&H
HC20K	-	X	-	-	-	KaVo
HC20KL	X	X	-	-	-	KaVo
HC21K	-	X	-	-	-	KaVo
HC21KL	X	X	-	-	-	KaVo
HC22K	-	-	X	-	-	KaVo
HC22KL	X	-	X	-	-	KaVo
HC8021K	-	X	-	-	-	KaVo

Article no.	Light	Standard head	Small head	Mini head	Direct connection	Quick-coupling system
HP21KL	X	X	-	-	-	KaVo
HP21NL	X	X	-	-	-	NSK
HP21SL	X	X	-	-	-	Sirona
HP21WL	X	X	-	-	-	W&H
HP22KL	X	-	X	-	-	KaVo
HP22NL	X	-	X	-	-	NSK
HP22SL	X	-	X	-	-	Sirona
HP22WL	X	-	X	-	-	W&H

KaVo, Sirona, W&H and NSK are registered trademarks. MK-dent does not have any economical connection to the companies mentioned above.

4.2 Combination with other products

REF MK-dent turbine	HB21K/KL, HB23K/KL, HE20K/KL, HE21K/KL, HE22K/KL, HC20K/KL, HC21K/KL, HC22K/KL, HC8021K, HP21KL, HP22KL	HE21SL, HE22SL, HP21SL, HP22SL	HE21W/WL, HE22W/WL, HP21WL, HP22WL	HE21N/NL, HE22N/NL, HP21NL, HP22NL
REF MK-dent coupling	QC4012K, QC4014K, QC5015K, QC5016K, QC5016KW, QC5116K, QC5116KW, QC6016K, QC6016KW, QC6116K, QC6116KW	QC5016SW, QC6016SW	QC4014W, QC5016W	QC4014NT, QC5016NT, QC6016NT
Connection	KaVo	Sirona	W&H	NSK

KaVo, Sirona, W&H and NSK are registered trademarks. MK-dent does not have any economical connection to the companies mentioned above.

4.3 Product components

Not applicable

4.4 Scope of delivery

The scope of delivery of the following components and accessories is included in the delivery unit:

HB12, HB14, HB15L, HB16L, HE11, HE12, HE14, HE15L, HE16L, HE17, HE17B, HE18			
AC0001	AC0003	LT1014P	--
Nozzle pin	Dental brush for maintenance of the spindle	Oil adapter for maintenance of the spindle	--

HB21K/KL, HB23K/KL, HE20K/KL, HE21K/KL, HE22K/KL, HC20K/KL, HC21K/KL, HC22K/KL, HC8021K, HP21KL, HP22KL			
AC0001	AC0003	LT1014P	LT1012P
Nozzle pin	Dental brush for maintenance of the spindle	Oil adapter for maintenance of the spindle	Oil adapter for maintenance of the rotor (KaVo)

HE21SL, HE22SL, HP21SL, HP22SL			
AC0001	AC0003	LT1014P	LT1017
Nozzle pin	Dental brush for maintenance of the spindle	Oil adapter for maintenance of the spindle	Oil adapter for maintenance of the rotor (Sirona)

HE21W/WL, HE22W/WL, HP21WL, HP22WL			
AC0001	AC0003	LT1014P	LT1018
Nozzle pin	Dental brush for maintenance of the spindle	Oil adapter for maintenance of the spindle	Oil adapter for maintenance of the rotor (W&H)

HE21N/NL, HE22N/NL, HP21NL, HP22NL			
AC0001	AC0003	LT1014P	LT1020
Nozzle pin	Dental brush for maintenance of the spindle	Oil adapter for maintenance of the spindle	Oil adapter for maintenance of the rotor (NSK)

When you open the delivery package for the first time, make sure that the seal on the packaging is undamaged and that all components are included in the scope of delivery. If anything is missing, contact your dealer immediately.

4.5 Required materials that are not included

NOTE

NOTE

The scope of delivery does not include a coupling.

The scope of delivery does not include an instrument holder, item no. RT2050.

The scope of delivery does not include a holding force measuring tool, item no. RT1020.

The scope of delivery does not include MK-dent Premium Service Oil LU1011.

4.6 Configurations

The configuration of the turbines is described in **4.1 Description** and **5.3 Spare parts** – only the connection options differ. The different connections are described as follows:

REF MK-dent turbine	HB21K/KL, HB23K/KL, HE20K/KL, HE21K/KL, HE22K/KL, HC20K/KL, HC21K/KL, HC22K/KL, HC8021K, HP21KL, HP22KL	HE21SL, HE22SL, HP21SL, HP22SL	HE21W/WL, HE22W/WL, HP21WL, HP22WL	HE21N/NL, HE22N/NL, HP21NL, HP22NL
Connection	KaVo	Sirona	W&H	NSK

KaVo, Sirona, W&H and NSK are registered trademarks. MK-dent does not have any economical connection to the companies mentioned above.

5 Accessories, consumables and spare parts

5.1 Accessories

Holding force measuring tool RT1020	
Instrument holder RT2050	
MK-dent Premium Service Oil LU1011	

5.2 Consumables

Not applicable

5.3 Spare parts

CAUTION

The use of unauthorized accessories or unauthorized modifications to the product can lead to injuries or instrument malfunctions.

→ Only use spare parts and components that have been approved by MK-dent for combination with the product.

	Turbine with direct connection
Eco Line replacement turbine	Standard head: HE11, HE12, HE15L, HE17, HE17B Small head: HE18 Mini head: HE14, HE16L
Basic Line replacement turbine	Standard head: HB12, HB15L Mini head: HB14, HB16L

	Turbine with coupling connection
Prime Line replacement turbine	Standard head: HC8021K, HP21KL/NL/SL/WL Small head: HP22KL/NL/SL/WL
Classic Line replacement turbine	Standard head: HC20K/KL, HC21K/KL Small head: HC22K/KL
Eco Line replacement turbine	Standard head: HE20K/KL, HE21K/N/W, HE21KL/NL/SL/WL Small head: HE22K/N/W, HE22KL/NL/SL/WL
Basic Line replacement turbine	Standard head: HB21K/KL Mini head: HB23K/KL

6 Installation

6.1 Requirements for the environment

To supply the products with clean and dry air, clean the compressor air filter in your dental practice. Do not supply moist compressed air. For drainage and regular maintenance of your compressor, please read the operating instructions for your air compressor. Ensure that the compressor is maintained in accordance with its instructions for use.

6.2 Unpacking

When opening the package for the first time, make sure that the content according to **4.4 Scope of delivery** is enclosed with the product. Make sure that the turbine has no defects or damage such as cracks.

6.3 Assembly

Not applicable

6.4 Preparation for use

WARNING

Danger of an infection for patients and users due to unsterile instruments
 → Before initial use and after each treatment, the instrument and accessories must be processed according to these instructions for use; see section **8 Processing according to ISO 17664-1**.

Malfunction and/or danger of infection from moist and contaminated compressed air and contaminated water for cooling
 → The product may only be used in the field of dentistry with a dental treatment unit that has a CE marking or is otherwise certified/approved according to the valid national legal regulations and that is provided with a suitable supply of water and compressed air in compliance with the applicable standards.

CAUTION

A damaged instrument or damaged components can injure the patient, the user and third parties.

- Test for proper functioning before each use and inspect the instrument visually for external damages.
- Do not use the instrument and components if there are any irregularities.

If a turbine is removed, water and oil will stay behind on the coupling connection. Remove the water and oil completely. Confirm that no water or oil remains in the coupling of your turbine hose. If there is still any residue left, remove it with a dry cloth.

Make sure that no water or oil remains in the connection of the instrument. Remove residues with a cotton swab.

6.4.1 Check and setting of the compressed air

WARNING

Excess or insufficient air pressure impairs the performance of the instrument and can cause pain in the patient. Problems can also arise during use, such as injuries due to instrument damage.

Connect a pressure gauge to the coupling on your turbine hose until you hear a locking sound. Make sure that the measuring instrument is securely attached. If you do not have a pressure gauge, contact your dealer.

Connect the instrument to the pressure gauge coupling until you hear a locking sound. Make sure that the measuring instrument is securely attached.

Check the following pressures (also see section **13 Technical data**):

- Drive air pressure
- Spray water pressure
- Spray air pressure
- Return air pressure

The compressed air must be set according to section **13 Technical data**. Remove the pressure gauge after setting the pressure.

6.4.2 Connecting to the coupling system

WARNING

Malfunction due to decoupling of the instrument during treatment.

An instrument that is not properly engaged can come loose from the hose coupling during treatment.


→ Check that the instrument is securely attached to the hose coupling by gently pulling on it before each treatment.

Damage due to inaccurate coupling.

Inaccurate coupling can destroy the high-pressure lamp or LED of the coupling or reduce its service life.

→ Check that the turbine is securely attached to the hose by gently pulling on it.

6.4.2.1 Connecting the instrument – turbine with direct connection

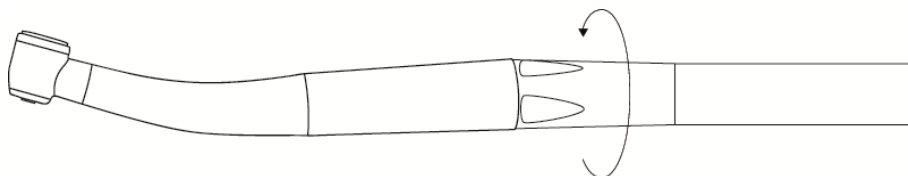- Place the instrument exactly onto the hose connection and turn the hose sleeve in the direction of the arrow until the instrument sits firmly.
- Check that the instrument is securely attached to the hose by gently pulling on it.
- For removal, hold the instrument firmly and unscrew the hose sleeve in the direction opposite to the arrow.

6.4.2.2 Connecting the instrument – turbine with coupling connection

- Place the instrument exactly onto the coupling and press firmly until you hear the connection snap in place.
- Check that the instrument is securely attached to the coupling by pulling on it.
- For removal, hold the instrument firmly and pull it off by turning it carefully.

6.5 Startup and operation

6.5.1 Before use

NOTE

NOTE

Before the start of every working day, the systems conducting water must be flushed for at least two minutes without any instruments attached. Another flushing process must be performed for 20 to 30 seconds after each patient.

6.5.2 Checking for correct operation

WARNING

Check the secure connection to the turbine hose after you hear a locking sound. An unstable connection can lead to injuries due to the instrument or other parts falling off.

Check for the absence of damage, cracks, breakage and corrosion on the product and that it is safe to use.

If you have problems, read section **11 Troubleshooting**. Use the water volume regulator on your treatment center to set the water supply.

6.5.3 Testing the water quantity

WARNING

Injury to the pulp and tooth.

A low quantity of spray water can result in overheating of the instrument and cause injuries.

→ Set the water quantity for the cooling water to at least 50 ml/min.

Check the water quantity with a measuring cup. Check the spray water tubes. If the water output is too low or a nozzle is clogged, clean the spray nozzle with the enclosed nozzle needle. Also see **8.1.2 Manual pre-cleaning of the spray nozzle**.

6.5.4 Using of rotating instruments

WARNING


- Risk of injury to patients if the rotating instrument falls out.
- After each use of rotating instruments, check the holding force of the collet chuck by pulling on the rotating instrument. Please be sure to do this while the drive is stationary.
- Risk of injury to patients and damage to the instrument if rotating instruments other than those listed below are used.
- Only use rotating instruments that comply with the data specified below.
 - Observe the instructions for use and proper use of the rotating instrument.
- Risk of injury due to the use of worn rotating instruments. Rotating instruments can fall out during treatment and injure the patient.
- Do not use rotating instruments with worn shafts.
- Risk of injury to users due to rotating instruments. Cuts and potential infections can be caused as a result.
- Use gloves or finger protection.

NOTE

NOTE

Use only rotating instruments made of steel or tungsten carbide that comply with DIN EN ISO 1797 (type 3) and meet the following criteria:

Diameter of the shaft:	1.590 to 1.600 mm
Full insert length of the shaft (standard head / small head / mini head):	min. 12/11/9 mm
Length of the shaft (standard head / small head / mini head):	max. 21/21/19 mm
Diameter of the cutter:	max. 2 mm

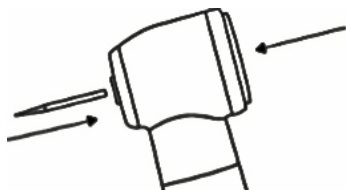- Press the button firmly with your thumb and simultaneously insert the rotating instrument as far as it will go.
- Pull on the drill to check that it sits firmly.
- Check the holding force of the collet at least once a week, using the RT1020 MK-dent dynamometer or another suitable device according to ISO 8325. The holding force of the collet must not fall below 2.5 kg. MK-dent recommends replacing the rotor when the holding force is below 3 kg.

7 Treatment

7.1 Patient, operator and assistant protection

WARNING

Before the first contact with the patient, the product must be properly processed to prevent infection. See section **8 Processing according to ISO 17664-1**.

Protect the patient appropriately before use.

During the treatment, the operator and assistants must wear gloves, safety goggles and masks to prevent infection.

CAUTION

Stop working whenever operating noises get louder or irregular, if there are stronger vibrations or in any case of overheating or visual damage. In such cases, contact your dealer or the MK-dent repair service.

MK-dent recommends lubricating the instrument after a treatment duration of 30 minutes, using MK-dent Premium Service Oil; see section **8.3.3 Maintenance, checks and testing**. If the instrument is to be used again directly after it is lubricated, it must run in idle for 15 - 20 seconds to remove excess oil. If the instrument is to be used again later on the same patient, retain it in an instrument holder (item no. RT2050) so that any excess oil can drain.


7.2 Removing rotating instruments

WARNING

Risk of injury from rotating instruments. Cuts and damage to the clamping system.

- Never touch rotating instruments.
- Never press the button while the instrument is rotating.
- Remove the rotating instrument after each treatment to prevent injuries and infections when putting it down.

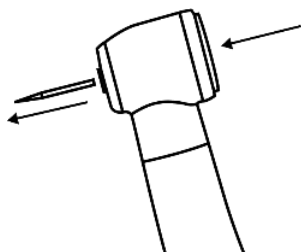- Use your thumb to press the button firmly after the rotating instrument has come to a complete stop and pull it out at the same time.

8 Processing according to ISO 17664-1

NOTE

NOTE

The lifecycle of the medical device is validated for 250 processing cycles. After 250 processing cycles, no significant negative effect on the product condition from processing was able to be detected. The possible number of applications does not depend on the processing cycles but rather on the application of the device and the wear and tear caused by this. The product's service life is therefore determined by signs of wear and damage occurring during use.

Dispose of damaged products in accordance with section **10 Transport and disposal**. Do not use any products with obvious signs of wear or damage.

Only use cleaning agents / disinfectants that are authorized for the processing of medical devices. Follow the manufacturer's instructions for the cleaning agents / disinfectants.

If unsuitable cleaning agents or disinfectants are used, this can have negative effects on the instruments:

- Damage or corrosion
- Discoloration of the product
- Corrosion of metal parts
- Reduced service life

Follow the instructions of the cleaning agent / disinfectant manufacturer. Do not use any abrasive cleaning agents / disinfectants.

The manufacturer recommends only using cleaning agents / disinfectants compliant with EN ISO 15883-1 and -2 for automated cleaning/disinfecting.

The instruments can be sterilized. Do not use any sterilization temperatures above 135 °C (275 °F).

Ultrasonic cleaning is not recommended for processing the instruments, as the ball bearings can be damaged, thereby reducing the product's service life.

8.1 Pre-treatment at the point of use

WARNING

Risk of infection from non-sterile instruments. Always wear protective gloves.

CAUTION

Using disinfecting baths and disinfectants containing chlorine can cause defects and malfunctions of the instrument.

The cleaning result is diminished by surface drying of residues.

→ Flush the transfer instrument immediately when you finish working on the patient

→ Immediately remove any cement, blood or composite residues

Do not submerge the instrument in solutions or similar.

NOTE

NOTE

Clean the coupling on your turbine hose and inside of the instrument connector according to **6.4 Preparation for use** before reconnecting the instrument.

- Flush the inner tubes of the product according to the RKI specifications for 20 seconds while it is still connected to the treatment unit.
- The instruments should be pre-cleaned directly after the treatment for the sake of personal safety. The aim is to prevent organic material and chemical residues from drying in the cavity or on the outer parts of the instruments and to avoid contaminating the environment.
- Remove adhering contaminants, body fluids and tissue using a disposable cloth / paper tissue.
- Flush the product with water after use (temperature below 35 °C (95 °F)).
- Do not use any abrasive cleaning agents or water with temperatures above 35 °C (95 °F), as this can lead to increasing adhesion of the contaminants, making the further processing steps more difficult.
- Clean the instrument within one hour after each treatment to prevent surface drying.
- Ensure that the instrument is stored safely and transported in dry condition in a closed container to the place of processing.

8.1.1 Preparation for cleaning

Take the drill out of the collet and pull off the instrument from the motor or quick coupling. Take off all removable parts (if any) and wash each part under running water.

Manual pre-cleaning is also recommended in the case of automatic processing, such as for removing blockages in the spray channels.


8.1.2 Manual pre-cleaning of the spray nozzle

WARNING

Injury to the pulp and tooth. A low quantity of spray water can result in overheating of the instrument and corresponding injuries.

→ Check spray water tubes and, if necessary, clean the spray nozzles with the nozzle pin (AC0001) included in the delivery.

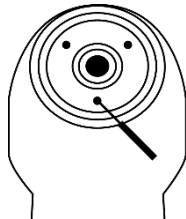Carefully remove clogging in the spray tubes with the enclosed nozzle pin.

8.1.3 Manual pre-cleaning of the instrument

CAUTION

Insufficient cleaning of the surface and water and air tubes.

Thorough manual pre-cleaning must be carried out to ensure successful cleaning.

The use of manual pre-cleaning requires training of the employees who carry out this processing instruction. This is necessary to ensure successful cleaning.

Tools:

- Drinking water¹ min. 100 ml 30°C ± 5°C (77 °F - 95 °F)
- Brush or soft toothbrush
- Dental brush (AC0003)
- 50 ml syringe or cannula
- Soft lint-free cloth

CAUTION


Risk of infection from dental brushes reused multiple times.

→ To avoid cross-infection, use a new dental brush for each instrument that is to be cleaned

Malfunctions can occur if cleaning is done in the ultrasonic machine.

→ Clean the instrument only as described below

¹ Drinking water quality according to the EU Drinking Water Directive (total germ count max. 100 CfU/ml)

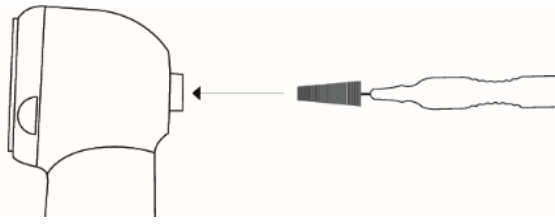- Carefully brush the complete instrument under running drinking water¹ using a brush or soft toothbrush to pre-clean the instrument.
- To pre-clean the collet, insert the enclosed dental brush (AC0003) into the collet and carefully loosen/remove any deposits and contaminants in rotating movements.
- Repeat this process several times and flush the dental brush with drinking water¹ each time.
- To pre-clean the water and air tubes, apply the syringe to the bottom end of the turbine and flush out the hollow spaces at least five times.

NOTE**NOTE**

Repeat these steps several times in case of heavy contamination.

8.1.4 Manual pre-cleaning of the glass rod of turbines with light

**CAUTION**

Risk of scratching and defects on the glass rod of the instrument if not handled with care.

→ Be extremely careful when you clean the glass rod

Conduct a visual inspection after each cleaning of the glass rod. Do not use the handpiece or angle handpiece if the glass rod is damaged and contact your dealer or the MK-dent repair service.

Wash the glass rod with warm drinking water¹ using a soft lint-free glass cleaning cloth.

Dry the glass rod using a compressed air gun or do so carefully using a soft lint-free glass cleaning cloth.

¹ Drinking water quality according to the EU Drinking Water Directive (total germ count max. 100 CfU/ml)

8.2 Manual processing

NOTE

NOTE

Use alcohol-based instrument disinfectant with a short residence time that does not cause any relevant adhesion of proteins on the instruments. The instrument disinfectant must be bactericidal, fungicidal and antiviral according to RKI requirements. Never use disinfectants containing chlorine, as these can damage the handpieces. Do not use aldehydic disinfectants, as they are incompatible with the material. Only use disinfectants whose microbiological efficacy is confirmed or guaranteed by the manufacturer (e.g., VAH/DGHM registration or CE marking).

According to the RKI recommendation, thermal disinfection in a steam sterilizer is required following the manual cleaning and disinfecting of medical devices of the “semi-critical B” class, especially if the instrument disinfectant used for the manual disinfection does not simultaneously have bactericidal, fungicidal and antiviral efficacy. Alternatively, the manual cleaning can serve as preparation for automatic cleaning.

Tools:

- ALPRO® WL-clean cleaning solution
- ALPRO® WL-cid disinfectant solution
- ALPRO® WL-dry drying spray
- ALPRO® cleaning adapter depending on the connector
- Lint-free disposable cloth

8.2.1 Manual cleaning

Connect the adapter to the Alpro® WL-clean spray bottle for cleaning the rotor and insert it into the turbine on the coupling side. Clean the instrument with 3 spray puffs for 2 seconds and with 1 minute resident time.

8.2.2 Manual disinfection, rinsing and drying

Disinfect the outer surfaces of the instrument using a disposable cloth with Alpro® WL-cid. Connect the adapter to the Alpro® WL-cid spray bottle for cleaning the rotor and insert it into the turbine on the coupling side and spray for 3 seconds. After 5 minutes of residence time, dry the instrument with Alpro® WL-dry spray.

CAUTION

The residence time after the last spray puff is 5 minutes.

8.3 Automatic processing

8.3.1 Automatic cleaning, disinfection and rinsing

NOTE

NOTE

MK-dent recommends thermal disinfection devices that satisfy the requirements of EN ISO 15883-1 and are approved by the manufacturer for dental handpieces and turbines. Cleaning should be performed at min. 55 °C (131 °F) for at least 6 minutes and disinfection at min. 90 °C (194 °F) for at least 5 minutes (for A0 value > 3000). A mildly alkaline cleaning agent with pH 9 to 11, such as Dr. Weigert Neodisher® MediClean forte, is recommended.

CAUTION

CAUTION

For automatic cleaning of the interior, use the adapters specified by the device manufacturer.

Only process instruments that are dry and free of residue using an automatic cleaning and disinfection process.

The automatic cleaning and disinfection process was validated using the device WD BHT INNOVA® M3 in program no. 03 (Thermal disinfection), with Dr. Weigert Neodisher® MediClean forte 0.5% as the cleaning agent. The instruments were connected with a silicon hose with a Luer-Lock adapter to the machine.

Parameters for automatic cleaning and disinfection:

- Pre-rinse for 1 minute with Neodisher® MediClean forte 0.5% cleaning agent at 30 °C (86 °F)
- Clean for 6 minutes with Neodisher® MediClean forte 0.5% cleaning agent at 55 °C (131 °F)
- Rinse for 1 minute with demineralized water
- Carry out thermal disinfection for 5 minutes at 90 °C (194 °F) using low-germ (< 10 CfU/ml) and low-endotoxin (< 0.25 endotoxin units/ml), demineralized water (for an A0 value > 3,000).

8.3.2 Automatic drying

The drying process is normally part of the cleaning process of your thermal disinfectant. Please follow the instructions for your thermal disinfectant. Allow the instrument to cool down before removing it from the adapter to ensure that this is safe and easier to do.

If there is still visible contamination in the thermal disinfectant, the process must be repeated.

To avoid any kind of damage to the instrument, make sure after each cycle that the instrument is dry inside and outside and free of residue. If moisture remains on or in the instrument after the cleaning process, dry it with a lint-free cloth. Also make sure that the cavities of the instruments are dry and dry them with medical compressed air if necessary.

8.3.3 Maintenance, checks and testing

CAUTION

Premature wear and malfunctions from improper maintenance and care
 → Reduced product service life
 → Perform appropriate maintenance and care at regular intervals


NOTE

NOTE

The instrument requires maintenance after each processing and before each sterilization. Use MK-dent Premium Service Oil LU1011 for this purpose.

The instrument must be maintained after each processing and before each sterilization. Carry out the following steps in the specified order:

- Screw the included oil adapter for maintenance of the spindle LT1014P onto the oil can. Spray into the spindle for 2 seconds.
- Screw the included oil adapter for maintenance of the rotor onto the oil can. Insert it into the instrument on the coupling side and spray oil into the drive air canal for 2 seconds.

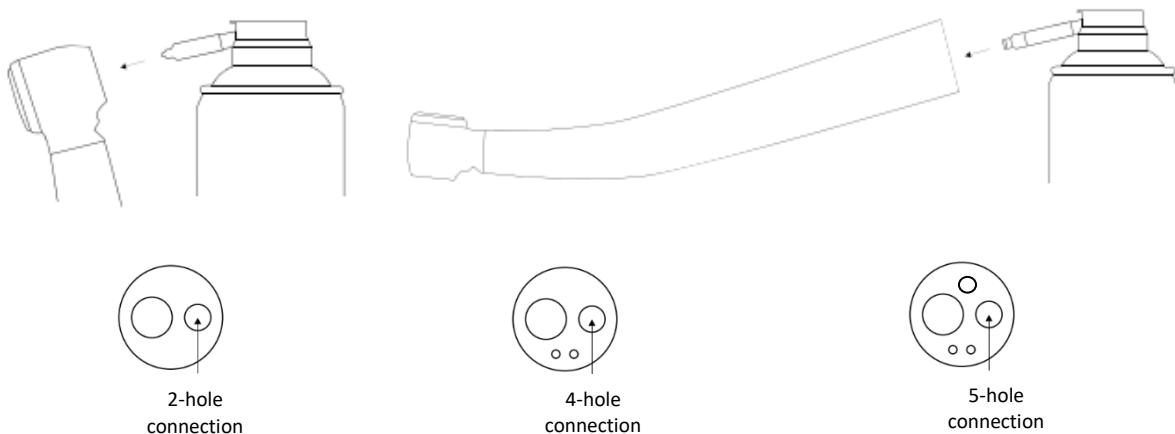8.4 Sterilization in an autoclave

8.4.1 Packaging

NOTE

NOTE

Use a validated packaging procedure that is suitable for sterilization. The sterilization pouch must be big enough for the instrument so that the packaging does not get stretched. The sterilization packaging must comply with EN ISO 11607-1 and be suitable for the sterilization process. Pack each instrument individually in a simple sterilization pouch. Seal the sterilization pouch carefully and check the seal for leaks.

8.4.2 Sterilization in a steam sterilizer according to EN 13060 / EN ISO 17665

WARNING

Premature wear and malfunctions from improper maintenance and care
 → Reduced product service life
 → Maintenance of the instrument with MK-dent Premium Service Oil is required before each sterilization.

Sterilization must be carried out before the first use and after each patient to prevent direct and cross-infection.

To avoid infection, wear suitable protective clothing such as a mask, gloves and safety goggles.

Contact corrosion from moisture can cause damage to the product.
 → Take the instrument out of the steam sterilizer directly after the sterilization cycle.

This medical device can withstand temperatures up to max. 135 °C (275 °F).

CAUTION

Depending on the conditions of use of your sterilizer or your sterilization methods, this product may undergo discoloration.

NOTE

NOTE

The sterilizer manufacturer's instructions for operation and load configuration must be followed explicitly. When sterilizing multiple instruments in one sterilization cycle, make sure that the instruments do not touch each other to ensure the best possible steam penetration. Observe the information in the sterilizer manufacturer's instructions for use.

After sealing it in the sterilization pouch, put the instrument on a sterilization tray to avoid direct contact with the inner chamber surface of the sterilizer and then place it into the autoclave. MK-dent recommends the use of a sterilizer in accordance with EN 13060, EN ISO 17665 or EN 285.

8.4.3 Sterilization parameters

Only use the following sterilization parameters to achieve an SAL of 10⁻⁶:

Autoclave with a fractioned pre-vacuum for at least 4 minutes at min. 134 °C (273 °F).

The sterilization was validated with the W&H LISA 517 device in a fractioned pre-vacuum process at 134 °C (273 °F) for 4 minutes. According to the sterilizer manufacturer's instructions for use, the drying time for the sterilization process is calculated by the device's software.

The instructions listed above for processing a medical device for reuse have been validated and found to be suitable by the medical device manufacturer. The user is responsible for ensuring that the processing actually performed, including the devices and materials used in the processing facility and any personnel involved, will achieve the intended result. This requires verification and/or validation and routine monitoring of the process.

8.4.4 After sterilization

When the sterilization cycle is finished, take the instrument out of the autoclave. When doing so, be sure to use appropriate personal protective equipment, as the instruments will be hot.

8.5 Storage

After sterilization and drying, check that there are no drops of residual water. Store the processed instruments in a dry, dark and cool room, protected from dust and with the lowest possible microbial load.

Temperature	-20 °C to +70 °C (-4 °F to +158 °F)
Relative humidity	5 % to 95 %, non-condensing
Air pressure	700 hPa - 1060 hPa (10 psi - 15 psi)

CAUTION

Instrument malfunction during operation after heavily cooled storage conditions
→ Allow heavily cooled instruments to warm up to room temperature before starting to use them.

The instrument must be completely dried before storage. Any remaining water can result in discoloration or rust.

Do not store the instrument in a place where it could be exposed to water.

Do not store the instrument in places where chemicals are stored or where corrosive gases exist.

9 Maintenance

9.1 Regular maintenance and checks

WARNING

Inappropriate maintenance or service can result in malfunctions and/or deterioration.

→ Always conduct proper maintenance and servicing according to section **8.3.3 Maintenance, checks and testing**

CAUTION

For safe use of the product, perform checks and maintenance of the instrument according to these instructions for use.

9.2 Cleaning in case of problems

In case of problems with the air or water supply lines, clean each part of the instrument. The instrument must be disconnected from the hose coupling for cleaning. Instructions can be found in section **8.1 Pre-treatment at the point of use**.

10 Transport and Disposal

CAUTION

Users or third parties may be infected by contaminated medical devices. Process the product before transport or disposal or dispose of it properly in compliance with the applicable regulations for the disposal of contaminated medical waste.

Process the instrument in accordance with section **8 Processing according to ISO 17664-1** and return it to us in the original packaging. We will take care of the professional disposal of the device and the original packaging for you.

11 Troubleshooting

NOTE

NOTE

This instrument must not be repaired by the user in the event of a malfunction. Please contact your dealer or the MK-dent repair service.

If this instrument does not function as described in these instructions for use, check the list below for your problem. If your problem is not resolved or if your problem is not listed below, please contact your dealer.

Symptoms	Possible causes	Countermeasures
1) Low efficiency	- No air supply - Compressed air is too low in the hose connection - The control system of the treatment unit is not working	- Clean according to 8.2.1 Manual cleaning - Set the motor pressure according to 13 Technical data - Test the control unit with a different handpiece - Request repair by your dealer
2) No water	- The water supply system of your treatment unit is switched off - No water gets into the turbine hose - Water tubes are clogged	- Activate the water supply system of your treatment unit - Check the water supply system of your treatment unit - Request repair by your dealer
3) Water leak near the turbine hose	O-ring(s) on your connection is (are) damaged.	- Replace the seals - Request repair by your dealer
4) No light	- Light guide is defective - Contactor of the light guide is defective	- Replace the glass rod - Repair/replace the connector - Request repair by your dealer
5) Irregular operating noises	- Ball bearing has run dry - Rotor imbalance - Ball bearing defective - Corrugated spring washers are worn out - O-rings are worn out	- Lubricate the turbine - Replace the rotor - Replace the corrugated spring washers - Replace O-rings - Request repair by your dealer
6) Strong vibrations	- Ball bearing defective - Ball bearing worn out	Request repair by your dealer

Symptoms	Possible causes	Countermeasures
7) Strong heat development	• Rotor has run dry • Insufficient water supply • O-rings are worn out • Ball bearing defective	• Lubricate/replace rotor • Clean/replace tube • Replace O-rings • Request repair by your dealer
8) Optical damages	• Drop damage • Deformation	• Replacement • Request repair by your dealer
9) Rotation speed inadequate	• Air supply insufficient	• Replace propellant air tube • Cleaning • Replace coupling • Request repair by your dealer

12 Repair

CAUTION

Contaminated medical devices can infect users or third parties. Process the instrument before dispatching it.

Process this instrument in accordance with section **8 Processing according to ISO 17664-1** and send it to your dealer or directly to MK-dent if it is defective. Please contact MK-dent customer service to arrange for shipping or collection for factory repair. Questions about noise development, maintenance, preparation and functionality can often be clarified in a conversation, making it unnecessary to send in the product.


13 Technical data

MK-dent recommends setting the spray air pressure at least 0.3 bar higher than the spray water pressure.

Article no.	Coating	Drive air pressure [bar]	Spray water pressure [bar]	Spray air pressure [bar]	Return air pressure [bar]	Air consumption [l/min]	Recommended pressure [N]	Spray water consumption [ml/min]	Sound level [dB(A)]	Idle rotation speed [min ⁻¹]
HB12	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	335,000
HB14	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	425,000
HB15L	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	335,000
HB16L	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	425,000
HB21K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	360,000
HB21KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	360,000
HB23K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	425,000
HB23KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 65	425,000
HE11	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	335,000
HE12	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	335,000
HE14	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	425,000
HE15L	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	335,000
HE16L	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	425,000
HE17	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	360,000
HE17B	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	360,000
HE18	Chrome	2.2-2.5	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	450,000
HE20K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	370,000
HE20KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	370,000
HE21K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	360,000
HE21KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	360,000
HE21N	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	360,000
HE21NL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	360,000
HE21SL	Chrome	2.7	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	360,000
HE21W	Chrome	3.0-3.3	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	370,000
HE21WL	Chrome	3.0-3.3	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	370,000
HE22K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	440,000
HE22KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	440,000
HE22N	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	440,000
HE22NL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	440,000
HE22SL	Chrome	3.0	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	450,000
HE22W	Chrome	3.0-3.3	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	475,000
HE22WL	Chrome	3.0-3.3	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 62	475,000
HC20K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 57	370,000
HC20KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 57	370,000
HC21K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 57	360,000
HC21KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 57	360,000
HC22K	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	440,000
HC22KL	Chrome	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	440,000
HC8021K	Titanium	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	360,000
HP21KL	Titanium	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 55	360,000
HP21NL	Titanium	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 55	360,000
HP21SL	Titanium	2.7	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 55	360,000

Article no.	Coating	Drive air pressure [bar]	Spray water pressure [bar]	Spray air pressure [bar]	Return air pressure [bar]	Air consumption [l/min]	Recommended pressure [N]	Spray water consumption [ml/min]	Sound level [dB(A)]	Idle rotation speed [min ⁻¹]
HP21WL	Titanium	3.0-3.3	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 55	370,000
HP22KL	Titanium	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 55	440,000
HP22NL	Titanium	2.8	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	440,000
HP22SL	Titanium	3.0	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	450,000
HP22WL	Titanium	3.0-3.3	0.8-2.5	1.0-4.0	< 0.5	40	2-3	> 50	< 60	475,000

14 Explanation of symbols

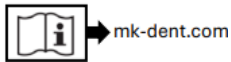	Caution – see section 1 Explanation of the warnings
	Sterilizable in a steam sterilizer (autoclave) up to the specified temperature
	Suitable for automatic processing in thermal disinfection devices
	Consult instructions for use and refer to the source of the electronic instructions for use
	CE marking with identification number of the Notified Body – indicates that the product is in compliance with European legislation for medical devices
	Medical device (CE conformity assessment procedure)
	Catalog number
	Serial number
	Rotor
	Manufacturer
	Date of manufacture
	Unique product identification – UDI (Unique Device Identifier) – indicates a barcode as containing Unique Device Identification
	HIBC Code (Health Industry Bar Code)

15 Guarantee and warranty

MK-dent assures you that your instrument has undergone a thorough **quality inspection** before leaving the factory. These quality inspections are carried out according to considerably stricter standards than the requirements of the product standard (ISO 14457 Dentistry - Handpieces and motors). In addition, we offer you comprehensive warranty services to ensure your satisfaction.

15.1 Guarantee

MK-dent grants the end customer a guarantee for perfect function and freedom from defects in material and workmanship from the date of purchase.

Exclusions:

- Wear parts, e.g., collets, ball bearings, gaskets, adapters, bulbs and O-rings
- Fall damage
- Willful damage and improper use contrary to the instructions for use
- Transport damage due to improper packaging (Instruments must be packed individually and with sufficient padding)
- No enclosure of the invoice copy of the device purchase (invoice recipient must be a dental practice)
- The warranty is valid from the date of purchase and is non-transferable
- Defects or their consequences due to interventions or modifications to the product by the end user or by a third party not authorized by MK-dent

Regions, time periods and validity:

International	
Basic Line	12 months
Eco Line	18 months
Classic Line	18 months
Prime Line	24 months

15.2 Legal warranty

MK-dent grants the statutory warranty on all products to the contractual partner. This covers material and processing defects that were already present at the time of the transfer of risk or that occur within the statutory warranty period. Warranty claims will only be recognized if the product is submitted with proof of purchase in the form of an invoice or a copy of the delivery note (contractual partner, date of purchase, model and serial number must be clearly visible).

Exclusions:

- Wear parts, e.g., collets, ball bearings, gaskets, adapters, bulbs and O-rings
- Fall damage
- Willful damage and improper use contrary to the instructions for use
- Transport damage due to improper packaging (instruments must be packed individually and with sufficient padding)
- The warranty is valid from the date of purchase and is non-transferable
- Defects or their consequences due to interventions or modifications to the product by the end user or by a third party not authorized by MK-dent

15.3 Conditions

- Proof of purchase (invoice or delivery note) is required for guarantee or warranty claims.
- The guarantee or warranty is void if the product is tampered with or modified by unauthorized third parties.
- Repairs or replacements under warranty do not entitle the warranty period to be extended or restarted.
- MK-dent reserves the right to repair or replace defective instruments.

15.4 Utilization

In the event of a warranty claim, please contact MK-dent customer service directly. In the event of a warranty claim, please contact your contractual partner.

Contact MK-dent customer service

Phone: +49 4532 40049-0

Email: service@mk-dent.com

Please have the following information ready:

- Product model and serial number
- Proof of purchase with date
- Detailed description of the problem

Our customer service team will guide you through the warranty process and ensure a quick and satisfactory solution. After the repair, your instrument will be returned to you securely packaged and sealed together with a detailed test report.

15.5 Exclusion of liability

Within the scope of the guarantee or warranty, MK-dent is not liable for financial losses, downtimes, loan or rental equipment, loss of profit or similar consequences. The liability of MK-dent is limited to the purchase value of the product.

MK-dent is not liable for any defects and their consequences that may arise due to natural wear and tear, improper handling, improper cleaning or maintenance, non-compliance with operating or handling instructions, calcination or corrosion, contaminated air and water supply or chemical or electrical influences that are unusual or not permitted according to the instructions for use of MK-dent or other manufacturer's instructions for use.

Any liability is excluded if defects or their consequences are caused by interventions or modifications to the product by the end user or by a third party not authorized by MK-dent.

MK-dent reserves the right to change the warranty conditions at any time. Changes have no retroactive validity for products already purchased.

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
5.1 Acessórios

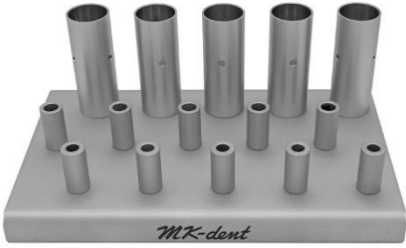Ferramenta de medição da força de retenção RT1020	 A digital force gauge with a black plastic housing and a stainless steel probe. The display shows '0.00'. Buttons for 'UNIT', 'TARE', 'LIGHT', and 'ON/OFF' are visible. The MK-dent logo is at the bottom of the housing.
Suporte de instrumentos RT2050	 A grey plastic base with 12 vertical slots of varying heights, designed to hold dental instruments. The MK-dent logo is printed on the front.
Premium Service Oil LU1011 da MK-dent	 A blue aerosol can with a white nozzle. The label features the MK-DENT logo, a drop icon, and the text 'PREMIUM SERVICE OIL'. Below this is a diagram of three dental turbines. At the bottom, it says 'LU1011 450', 'Volume = 500ml (l)', and 'NET 12.25 oz (349g)'.

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